☆☆☆
TRISTAR

🍷 🌊 🐾 🧑
WG-2424

FIT

EN **INSTRUCTION MANUAL**

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PT **MANUAL DE UTILIZADOR**

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SV **BRUKSANVISNING**

TRISTAR

Dear customer,

Congratulations and thank you for purchasing this high quality product. Please read the instruction manual carefully so that you will be able to make the best possible use of the appliance. This manual includes all the necessary instructions and advice for using, cleaning and maintaining the appliance. If you follow these instructions you will be guaranteed of an excellent result, it saves you time and will avoid problems. We hope you will have a great deal of pleasure from using this appliance.

**YOU CAN FIND ALL INFORMATION AND SPARE PARTS
WWW.SERVICE.TRISTAR.EU!**

	An answer to all your questions
	Tips and tricks on how to use your product
	Order spare parts for your product online
	Register for software updates for your product

SAFETY INSTRUCTIONS

- Please read the instruction manual carefully before using the device. Please keep these instructions, the guarantee certificate, the sales receipt and, if possible, the carton with the inner packaging.
- This appliance is not intended to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instructions concerning the use of the appliance by a person that is responsible for their safety.

- By ignoring the safety instructions the manufacturer can not be hold responsible for the damage.
- To protect children against the dangers of electrical appliances, please make sure that you never leave the appliance unattended. Therefore you have to select a storage place for the appliance where children are not able to grab it. Make sure that the cable is not hanging in a downward position.
- This appliance is only to be used for household purposes and only for the purpose it is made for.
- The appliance must be placed on a stable, level surface.
- Handle the balance with care, as it is a precise measuring instrument. Do not drop it or jump on it.
- Do not load the appliance with more than 180kg to avoid damaging the device.
- If the scales are carried or not going to be used, please store them only in a horizontal position.

THE BODY FAT ANALYSIS SCALE IS NOT SUITABLE FOR THE FOLLOWING PERSONS:

- Children under 10 years and adults over 85 years.
- Persons with symptoms of fever, oedema, or osteoporosis.
- Persons undergoing treatment for dialysis.
- Persons who use a pacemaker.
- Persons who take cardiovascular medication.
- Pregnant women.
- Sportspersons who perform over 5 hours of intensive sporting activities each day.
- Persons with a heartbeat lower than 60 beats/minute.

BATTERIES

- Do not expose the battery to thigh temperatures or direct sunlight. Never throw batteries into the fire. There is a danger of explosion!
- Keep batteries away from children. They are not a toy!

TRISTAR EN INSTRUCTION MANUAL


- Do not open the batteries by force.
- Avoid contact with metallic objects.(Rings,nails,screws et cetera.)there is a danger of short-circuiting!
- As a result of a short-circuit batteries may heat up consider-ably or even catch fire.this may result in burns.
- For your safety the battery poles should be covered with adhesive strips during transport.
- If a battery leaks do not rub the liquid into contact with the liquid they should be washed thoroughly and eyes rinsed out with clear water,if symptoms persist please consult a doctor."

PARTS DESCRIPTION

1. Contact surface
2. Display
3. Down button
4. Up button
5. Set button

DISPLAY SYMBOL	MEANING
Kg/lb/st)	Body weight
P0-9	User memory space
	Male
	Female
CM(FT;In)	Body height
Age	Age
FAT	Body fat display
TBW	Body water display
MUS	Muscle mass diaplay
BONE	Bone weight display
%	Percentage-body fat/body water/muscle mass/ bone weight

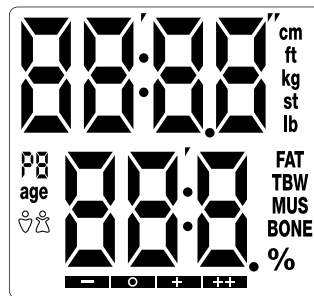BEFORE THE FIRST USE

- Take the appliance and accessories out the box. Remove the stickers, protective foil or plastic from the device.
- The device can only function correctly if you stand on the contact surface with bare, dry feet!
- The main principle of this device is based on the measurement of the electric impedance within the human body. In order to determine the actual body data, the device leads a weak, unnoticeable electric signal through your body.
- This signal is generally safe! However, please read our special safety instructions!
- The method of measurement is referred to as "bioelectric impedance analysis" (BIA), a factor which is based on the relationship between body weight and associated body water percentages as well as on other biological data (age, gender, body weight).
- There are heavy persons with a well-balanced percentage between fat and muscle mass, and there are seemingly slim persons with significantly high body fat values, fat is a factor that puts serious strain on the body, in particular the cardio-vascular system. Therefore, controlling both values (fat and weight) is a prerequisite for health and fitness.

USE

- The scales are factory set with the units "kg" and "cm". To change the settings to "lb" or "kg" and "Ft:in", proceed as follows: On the back side of the scale beside the battery compartment a switch can be found, push the below switch, choose between the measuring units kg/lb/st.

WEIGHING FUNCTION

- Place the scale on a hard, flat surface (avoid carpets). An uneven floor contributes to inaccurate measurements.
- Stand on the scales with both feet. The device switches on automatically.
- Distribute your weight equally and wait for the scale to determine your weight.
- The digits on the display will flash before your exact weight will be shown.
- In order to extend the battery life, the device will automatically turn off 10 seconds after you have stepped off the scale.

SAVING PERSONAL DATE

- The scale has default values that you can change, as described below.
- Carry out the following steps in sequence and within a few seconds, otherwise the device will automatically cancel the input mode.
- Press the button SET. The memory space will flash on the display (P0). Confirm the setting by pressing SET or select another memory space by pressing the UP/DOWN button and then SET to go to the next setting.
- The symbol for "male" or "female" will flash. Confirm the setting by pressing SET or use the UP/DOWN button to choose the symbol for your gender and then press to go the next setting.
- The height setting will flash. Select your height by pressing the UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.
- The setting AGE will flash. Select your age by pressing the UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.
- The device can only function correctly if you stand on the device with bare, dry feet. If necessary wipe the contact surface to improve the results.
- Due to the significant fluctuation of body fat the optimal testing time is between 6 p.m. and 8 p.m., since the body fat remains very consistent during this period of time.

- Press the SET button, the user number will flash (P0-P9). Use the UP/DOWN button to select your preset data.
- Wait for the user number to stop flashing.
- Stand on the contact surfaces with both feet.
- Firstly your weight will flash on the display.
- After a few seconds, the following results will be shown:
 - the body fat content in% (FAT symbol)
 - the results analysis in the form of a bar
 - the body water content in % (TBW symbol)
 - the muscle mass content in % (MUS symbol)
 - the bone mass content in %(BON symbol)

ERROR

If an error message appears on the display, wait until the device is turned off. Turn it on again to restart the operation. If the error display has not disappeared, take the batteries out of the battery compartment for 10 seconds.

- 8888 = Self test (when switching on) to check if the numeric display functions properly.
- `---- = The scale is not placed properly.
- Err = The percentage of body fat content lies above or below the limits of the measuring scale
- EEEE = The weight lies above the allowed maximum weight.

CLEANING AND MAINTENANCE

- Do not immerse the appliance in water or any other liquid.
- Clean the appliance with a damp cloth. Never use harsh and abrasive cleaners, scouring pad or steel wool, which damages the device.

GUARANTEE

- Tristar is not liable for damages occurred through:
 - o In case the appliance has fallen.
 - o In case the appliance has been technically changed by the owner or another third party.
 - o In case of improper use of the appliance.
 - o In case of normal wear of the appliance.
- By executing repairs the original warranty period of 2 years will not be extended, nor the right to a complied new warranty. This warranty is only legal on European soil. This warranty does not overrule the European directive 1944/44CE.
- Always keep your receipt, without this receipt you can not claim any form of warranty whatsoever.
- Damage caused by not following the instruction manual will lead to a void of warranty, if this results in consequential damages Tristar will not be liable.
- Tristar will not be liable for material damage or personal injury caused by improper use or if the safety instructions are not properly executed.
- Cleaning, as mentioned in this manual, is the only necessary maintenance for this appliance.
- When the appliance should be repaired, make sure that this will be performed by an authorized company.
- This appliance may not be amended or changed.
- If problems arise during the 2 years from the date of purchase, which are covered by the factory warranty, you can go back to the point of purchase to replace it for a new one.
- For questions or compliances please contact your retailer "point of purchase."
- This appliance is covered with a 24 month guarantee starting on the date of purchase (receipt).

- Only material or manufacturing defects are included in this warranty.
- If you wish to make a claim please return the entire machine in the original packaging to your dealer together with the receipt.
- Damage to accessories does not mean automatic free replacement of the whole appliance. In such cases please contact our hotline. Broken glass or breakage of plastic parts are always subject to a charge.
- Defects to consumables or parts subject to wearing, as well as cleaning, maintenance or the replacement of these parts are not covered by the guarantee and hence are to be paid!
- The guarantee lapses in case of unauthorized tampering.
- After the expiry of the guarantee repairs can be carried by the competent dealer or repair service against the payment of the ensuing costs.

GUIDELINES FOR PROTECTION OF THE ENVIRONMENT

This appliance should not be put into the domestic garbage at the end of its durability, but must be offered at a central point for the recycling of electric and electronic domestic appliances. This symbol on the appliance, instruction manual and packaging puts your attention to this important issue. The materials used in this appliance can be recycled. By recycling used domestic appliances you contribute an important push to the protection of our environment. Ask your local authorities for information regarding the point of recollection.

PACKAGING

The packaging is 100% recyclable, return the packaging separated.

PRODUCT

This device is equipped with a mark according to the European Directive 2012/19/EU on Waste Electrical and Electronic Equipment (WEEE). By ensuring that this product will be correctly processed into waste, it will help prevent possible negative consequences for the environment and human health.

EC DECLARATION OF CONFORMITY

This device is designed, manufactured and marketed in accordance with the safety objectives of the Low Voltage Directive "No 2006/95/EC", the protection requirements of the EMC Directive 2004/108/EC "Electromagnetic Compatibility" and the requirement of Directive 93/68/EEC.

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