MELISSA

16770048

DK	Massagepistol	2
NO	Massasjepistol	9
SE	Massagepistol	16
FI	Hierontavasara.....	23
UK	Massage gun	29
DE	Massagepistole	36

www.adexi.eu

INDLEDNING

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INTRODUCTION

To get the best out of your new massage gun, please read these instructions carefully before using it for the first time.

Take particular note of the safety precautions. We also recommend that you keep the instructions for future reference, so that you can refer to them at a later date.

SAFETY INSTRUCTIONS

- This appliance may only be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they are supervised or have been given instructions concerning use of the appliance in a safe way and understand the hazards involved.

- Children must not play with the appliance. Cleaning and user maintenance must not be performed by children without supervision.
- Never leave unattended when in use. Keep an eye on any children in the vicinity of the appliance when it is in use. The appliance is not a toy.
- Incorrect use of this appliance may cause personal injury or damage the appliance.
- Use for the intended purpose only. The manufacturer is not responsible for any injury or damage resulting from incorrect use or handling (see also Warranty Terms).
- For domestic use only. Not for outdoor or commercial use.

UK

- Remove all packaging and transport materials from the inside and outside of the appliance.
- Check that the appliance has no visible damage and that no parts are missing.
- Check that it is not possible to pull or trip over the cord of the appliance's USB cable during charging.
- Turn off the appliance before cleaning and when not in use.
- Keep the USB cable and appliance away from heat sources, hot objects and naked flames.
- The USB cable must not be twisted or wound around the appliance.
- Do not use the appliance if it has been dropped on the floor or in water, or has been damaged in any other way.
- The appliance may only be connected to a 5 V USB power supply using the USB cable supplied. The warranty is not valid if it is connected to the wrong voltage.
- Never try to repair the appliance yourself. Please contact the store where you bought the appliance for repairs under warranty.
- Unauthorised repairs or modifications will invalidate the warranty.
- Do not use accessories or spare parts which do not come from the supplier.
- Not for use in the bath.

UK

- Never place the appliance anywhere it can fall or be pulled into water, e.g. a wash basin or bath tub.
- Do not use the appliance for long periods at a time, as this can cause overheating and shorten the appliance's lifetime. In the event of overheating, switch the appliance off immediately and allow it to cool down completely before using it again.

MEDICAL WARNINGS

- **The appliance is not designed for physiotherapeutic treatment, but only for 'well-being' massage, relaxation of tense muscles and/or relief of low-level muscular pain.**

Use of the appliance cannot substitute professional physiotherapeutic or medical treatment.


- **If unsure about your medical condition, consult your doctor before using the appliance.**
- **The appliance must not be used by anyone with a pacemaker.**
- **The appliance must not be used by anyone with heart problems or pains of unknown origin or if pregnant.**
- **Do not use the appliance if you have taken painkillers, relaxants or alcohol, as your tolerance of heat can be weakened.**
- **Do not use the appliance from the neck upwards, and on sensitive areas.**

- **Do not use the appliance directly on swollen, inflamed or damaged skin areas.**
- **Do not use the appliance together with oils, creams or similar.**
- **Do not use the appliance just before going to bed, as the appliance has a stimulating effect that can make it difficult to fall asleep.**
- **Do not use the appliance on the same part of the body for too long; move the appliance over a larger area of the body in a continuous movement.**

UK

MAIN COMPONENTS


1. Handle
2. Socket
3. Massage head
4. Control panel
5. On-off / Gear increase button
6. Gear display
7. Gear decrease button
8. Charging socket
9. Charging indicator light

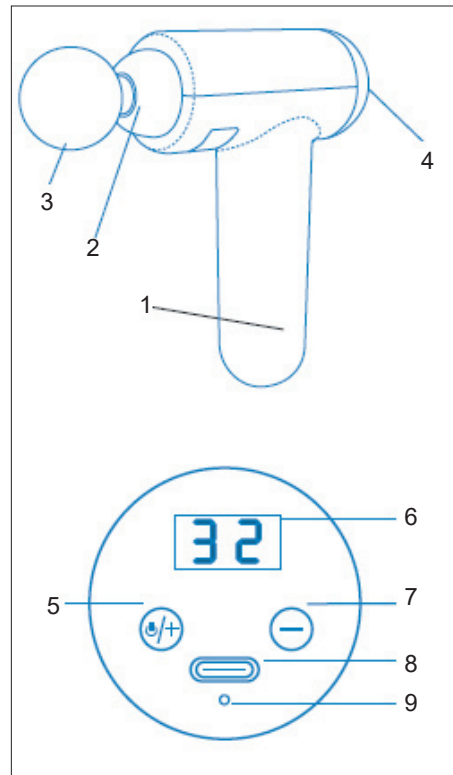USE


Massage heads

- A. Round massage head for large muscle areas, hands, back, buttocks and legs
- B. Spherical massage head for deep tissue massage and massage of joints, palms and soles
- C. Flat massage head for all muscle areas and muscle relaxation
- D. Fork-shaped massage head for Achilles tendons and muscles of the neck and spine

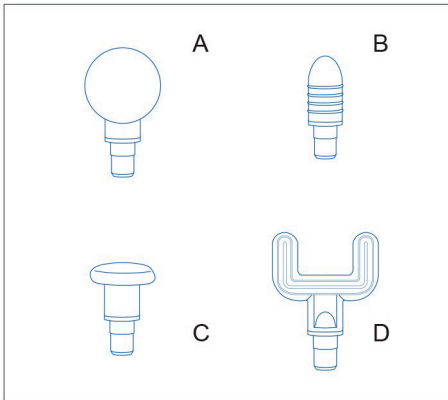The appliance helps to increase blood circulation, reduces muscle soreness/ aches and promotes recovery after exercise.

1. Insert the massage head (A/B/C/D) into the socket (2) to attach it onto the appliance.

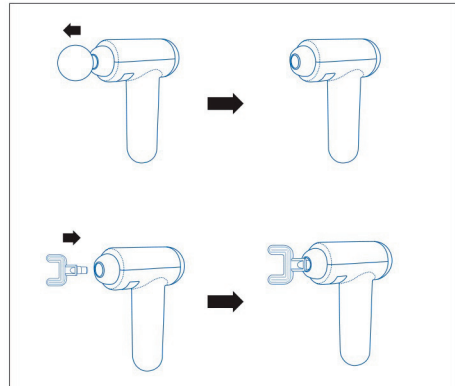2. Switch on the appliance by long pressing the on/off button (5). The appliance starts on default at the gear 10 setting. The number of the gear setting is displayed.
3. Press the increase button (5), or the decrease button (7) repeatedly to adjust the gear setting to have the massage intensity desired. There are 32 gears totally.
4. Press the massage head (3) against the skin with a firm grip on the appliance handle (1).
5. Each massage process will automatically stop in 15 minutes. Or stop the massage by long pressing on the On/off button (5).

CHARGING

Use the USB cable provided to connect the appliance to a 5V DC USB power output to charge the appliance.

A full charge takes approximately 3 hours and provides 60 minutes of operating power.

The indicator light (9) comes on red during charging and will turn green after completion.

NOTE: THIS PRODUCT CAN NOT BE CHARGED BY FAST CHARGING ADAPTOR. DON'T USE THE APPLIANCE DURING CHARGING.

NOTE: ALWAYS LET THE APPLIANCE COOLING DOWN BETWEEN 2 MASSAGE CYCLES. NEVER RUN THE APPLIANCE FOR TOO LONG TIME.

CLEANING

- The appliance must never be exposed to or immersed in water.
- Do not use scouring powder, steel wool or other aggressive cleaning materials to clean the appliance, as they can scratch the surface of the device.
- Use a soft, damp cloth to wipe the appliance instead.
- The removable massage heads can be washed in warm soapy water.
- Make sure they are completely dry before fitting them on the appliance again.

UK

INFORMATION ON DISPOSAL AND RECYCLING OF THIS PRODUCT

Please note that this Adexi product is marked with this symbol:

This means that this product must not be disposed of along with ordinary household waste, as electrical and electronic waste must be disposed of separately.

According to the WEEE directive, every member state must ensure correct collection, recovery, handling and recycling of electrical and electronic waste. Private households in the EU can take used equipment to special recycling stations free of charge.

In some member states you can in certain cases return used equipment to the retailer from whom you purchased it, if you are purchasing new equipment. Contact your retailer, distributor or the municipal authorities for further information on what to do with electrical and electronic waste.

WARRANTY TERMS

The warranty does not apply:

- if the above instructions have not been followed
- if the appliance has been interfered with
- if the appliance has been mishandled, subjected to rough treatment, or has suffered any other form of damage
- if the appliance is faulty due to faults in the electricity supply.

Due to the constant development of our products in terms of function and design, we reserve the right to make changes to the product without prior warning.

FREQUENTLY ASKED QUESTIONS

If you have any questions regarding the use of the appliance and cannot find the answer in this user guide, please try our website at www.adexi.eu.

You can also see contact details on our website if you need to contact us for technical questions, repairs, accessories or spare parts.

MANUFACTURED IN CHINA FOR

Adexi A/S

Lægårdsvej 9C

DK-8520 Lystrup

Denmark

www.adexi.eu

We cannot be held responsible for any printing errors

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